




Aerodynamics

AF27

Nozzle Flow Apparatus

Demonstrates the thermodynamics and fluid mechanics of the adiabatic expansion of air through subsonic and supersonic nozzles

Works with
VDAS®

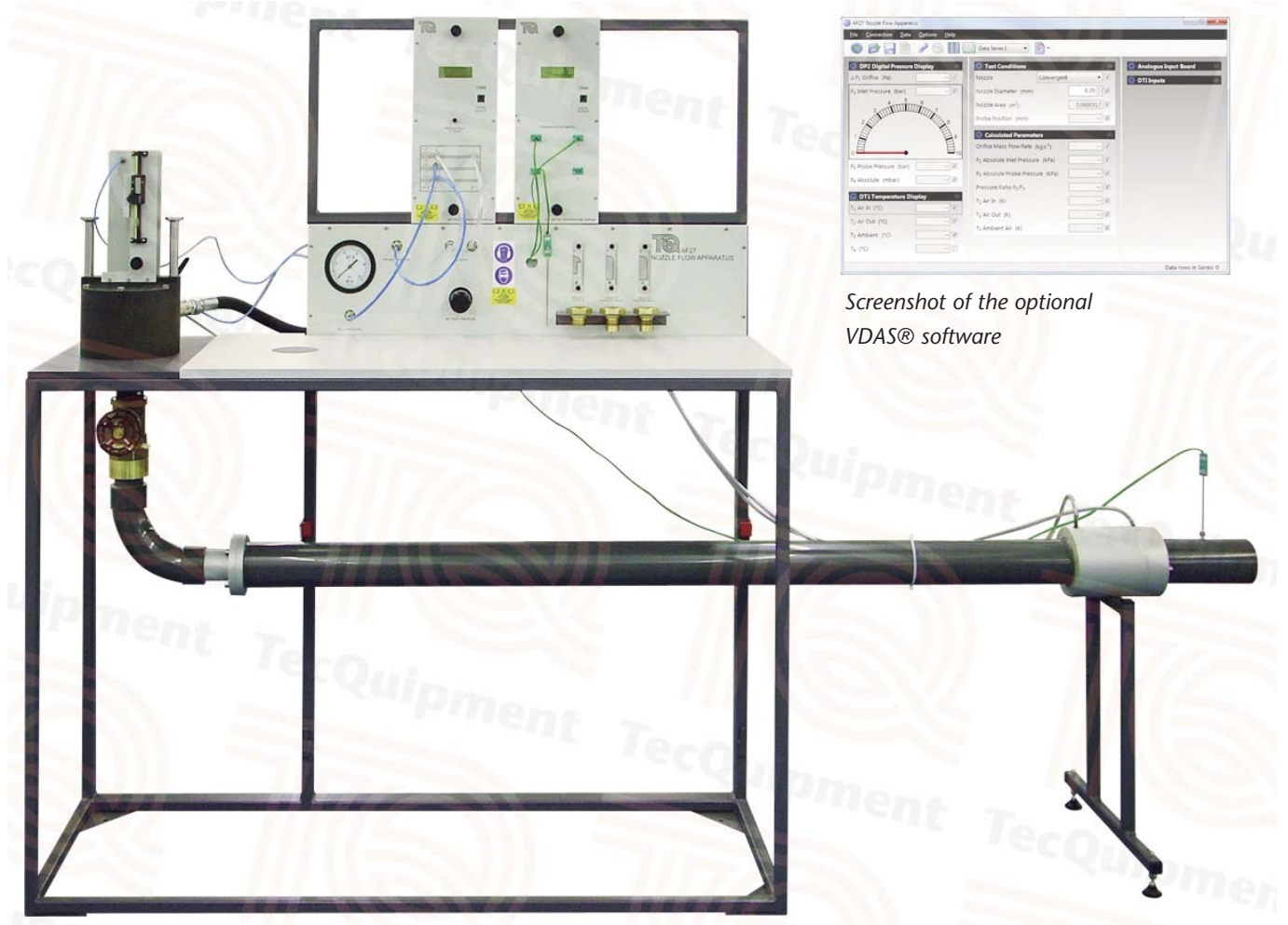

Screenshot of the optional VDAS® software

- Floor-standing apparatus that demonstrates the thermodynamics and fluid mechanics of the adiabatic expansion of air through subsonic and supersonic nozzles
- Connects to suitable laboratory compressed air supply or TecEquipment's optional Compressor (AF27a)
- Includes three interchangeable, profiled and polished brass nozzles – convergent, convergent-divergent and convergent-parallel
- Built-in instrument frame and worktop
- Electronic instruments measure and display multiple pressures and temperatures at the same time, for ease of use and for connection to TecEquipment's VDAS®
- Works with TecEquipment's Versatile Data Acquisition System (VDAS®) for instant recording of multiple readings and automatic calculations

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecquipment.com • **W** www.tecquipment.com
- An ISO 9001 certified company
- VDAS is a registered trademark of TecEquipment Ltd

TecEquipment

AF27

Nozzle Flow Apparatus

Description

The apparatus connects to TecEquipment's optional Compressor (AF27a) or a suitable laboratory supply of dry, clean compressed air.

It demonstrates the thermodynamic and fluid properties of the adiabatic expansion of subsonic and supersonic air flow through nozzles.

Its floor-standing frame holds:

- a pressure chest with a removable lid and nozzle traverse mechanism;
- a useful worktop;
- a pressure regulator to maintain the inlet/upstream pressure, with an analogue reference pressure gauge;
- three interchangeable, profiled and polished brass nozzles with mimics that fit on the traverse mechanism;
- an instrument frame with digital pressure and temperature displays

Students fit a nozzle into the chest (you test one nozzle at a time). Compressed air passes through the pressure regulator and an isolating valve. It then enters the pressure chest and passes vertically down through the nozzle, then through a precision downstream valve. The airflow then settles as it passes along a horizontal pipe, through an orifice and out to atmosphere.

The temperature and pressure displays accurately measure temperatures and pressures at key points around the apparatus, including the pressures around the orifice which students use to determine overall mass flow.

For ease of visibility and for good engineering practice, the analogue gauge also shows pressure in the chest.

A stainless-steel probe on a manually adjustable, vertical traverse measures the pressure distribution along the axis of the nozzle. A digital indicator measures the probe position in the nozzle. The traverse mechanism includes a pointer and a mounting for a 'mimic' of each nozzle. The pointer moves along the mimic to help students visualize the position of the probe.

The instrument frame has extra space for the optional VDAS® interface unit.

VDAS® allows accurate real-time data capture, monitoring, display, calculation and charting of all the important readings on a suitable computer (computer not included).

Standard Features

- Supplied with comprehensive user guide
- Five-year warranty
- Made in accordance with the latest European Union directives

Essential Ancillaries

- Compressor (AF27a)
- Versatile Data Acquisition System – Frame-mounted version (VDAS-F)

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecequipment.com • **W** www.tecequipment.com
- An ISO 9001 certified company
- VDAS is a registered trademark of TecEquipment Ltd

Experiments

- The relationship between pressure ratio and flow for convergent and convergent/divergent Laval nozzles
- The pressure profile in convergent/divergent nozzles at various pressure ratios
- Investigation of expansion with friction in a parallel passage at high subsonic velocities
- Boundary layer growth under subsonic and supersonic conditions
- The phenomenon of choked flow corresponding to sonic velocity at a nozzle throat

Essential Services

Electrical supply:

50 Hz to 60 Hz 90 VAC to 250 VAC for the instruments of the Nozzle Flow Apparatus.

Supply for the optional Compressor (AF27a) is determined on order

Air supply needed:

Clean, dry air at 2.5 kg.min⁻¹

Minimum 7 bar, Maximum 10 bar

Air connection: 3/4" BSP

Floor space needed:

3.9 m x 2 m for the Nozzle Flow Apparatus (including 1 m of free space in front of outlet)

Nominally 2.5 m x 1.5 m for the optional Compressor (AF27a)

Operating Conditions

Operating environment:

Laboratory

Storage temperature range:

-25°C to +55°C (when packed for transport)

Operating temperature range:

+5°C to +40°C

Operating relative humidity range:

30% to 95% (non-condensing)

Sound Levels

This equipment emits sound levels greater than 100 dB(A). You must wear ear defenders when you use it or work near to it.

Specifications

Dimensions and weight:

Nett: Length 2900 mm x depth 600 mm x height

1730 mm and 123 kg (including instruments)

Packed: 4.0 m³ and 218 kg



TecEquipment

tradition.

innovation.

integration.

infoWERK is a leading expert in the development of eLearning courseware, learning system solutions, teaching and AV equipment.

Furthermore infoWERK is the representative and system integrator of "TecQuipment".

TecQuipment is one of the global leaders in technical teaching equipment for engineering. If you are interested in one of TecQuipment's products feel free to contact us at:



infoWERK Medien & Technik GmbH

Martinsbühel 6 / A-6170 Zirl / Austria

Phone: +43 (0) 5238 52099-0 / Fax: +43 (0) 5238 52099-40

E-Mail: info@infowerk.at / Website: infowerk.at

Otto-Dürr-Straße 25

D-70435 Stuttgart, Zuffenhausen/ Germany

Phone: +49 (0) 711 342471-0 / Fax: +49 (0) 711 342471-11

E-Mail: info@de.infowerk.at / Website: infowerk.at