


Aerodynamics

AF300

Intermittent Supersonic Wind Tunnel

Investigates subsonic and supersonic air flow, including flow around two-dimensional models

Works with
VDAS®

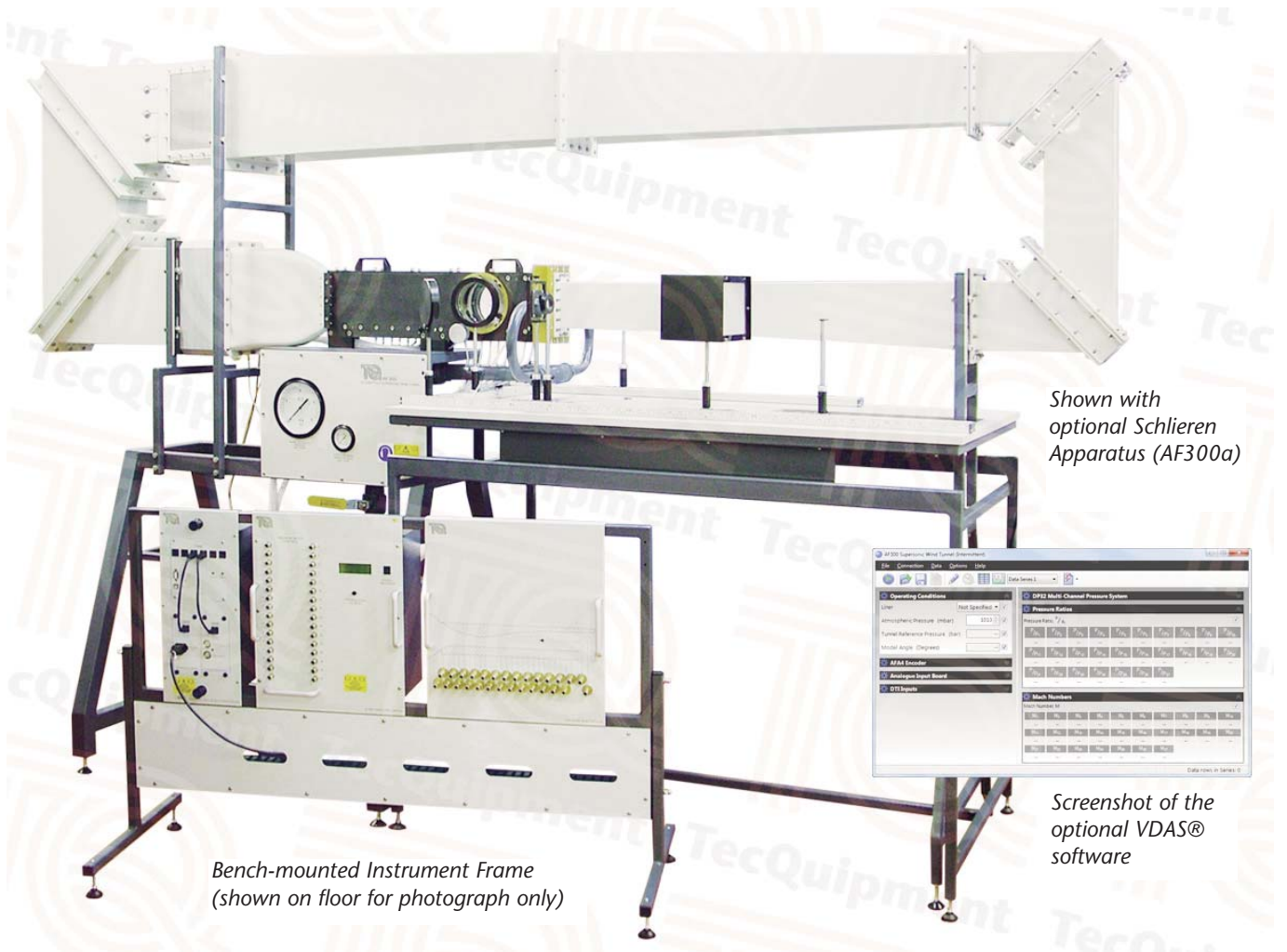

Shown with optional Schlieren Apparatus (AF300a)

Bench-mounted Instrument Frame (shown on floor for photograph only)

Screenshot of the optional VDAS® software

- Laboratory-scale wind tunnel for subsonic and supersonic tests, up to Mach 1.8
- Compact design – does not need large laboratory space
- Supplied with aerodynamic models for supersonic tests – includes model angle-feedback encoder
- Works with TecEquipment's Versatile Data Acquisition System (VDAS®) for instant recording of multiple readings and automatic calculations
- Electronic instruments measure and display multiple pressures at the same time, for ease of use and for connection to TecEquipment's VDAS®
- Supplied with set of different liners for controlled subsonic and supersonic air flow
- Induction flow for better air flow and accurate results

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecquipment.com • **W** www.tecquipment.com
- An ISO 9001 certified company
- VDAS is a registered trademark of TecEquipment Ltd

AF300

Intermittent Supersonic Wind Tunnel

Description

An intermittent operation, induction-type supersonic wind tunnel for investigations into subsonic and supersonic flow. This includes tests on the flow around two-dimensional models at subsonic and supersonic air speeds.

A compressed air supply (AF300b, available separately) induces a flow in the working section of the wind tunnel. This gives a less turbulent and more stable flow for accurate results and comparison with theory. The optional compressed air supply includes filters and air dryers to give a dust-free and dry air source needed for good results.

Students use a delivery valve to allow compressed air to enter the wind tunnel. The wind tunnel includes two analogue pressure gauges. One measures the compressed air pressure available from the supply (for reference); the other measures the pressure delivered to the wind tunnel and includes an electronic transducer that connects to TecEquipment's optional Versatile Data Acquisition System (VDAS®) to record the pressure.

The working section of the wind tunnel is a convergent-divergent nozzle with a removable top part ('liner'). The shape of the liner controls the maximum air velocity at the divergent part of the working section. Included are three different liners.

High optical-quality glass windows ('portals') are at each side of the divergent part of the working section. The portals allow students to use the optional Schlieren Apparatus (AF300a, available separately). This allows display and recording of images of pressure waves around two-dimensional models.

Included is a set of two-dimensional models. These mount between the portals inside the working section. Students can adjust the angle of the models. An encoder electronically measures the model angle.

Spaced at precise intervals along the working section of the wind tunnel are pressure tapings. Two extra tapings connect to one of the models when in use. A 32-way pressure display (included) connects to all the pressure tapings. It displays the pressures and transmits them to VDAS® for instant recording and calculations of pressure ratios and Mach numbers.

Included is a bench-mounting instrument frame that holds and provides power for the electronic instruments and the optional VDAS® interface unit. The instrument frame connects to a suitable electrical supply.

VDAS® allows accurate real-time data capture, monitoring, display, calculation and charting of all the important readings on a suitable computer (computer not included).

Standard Features

- Supplied with comprehensive user guide
- Five-year warranty
- Manufactured in accordance with the latest European Union directives

Experiments

- Pressure distribution along a convergent/divergent (Laval) nozzle with subsonic and supersonic air flow
- Comparison of theoretical and actual pressure distribution
- Comparison of actual and theoretical area ratio of a nozzle at supersonic air velocities (Mach numbers)
- Pressures around a two-dimensional model in subsonic and supersonic flow conditions, at different angles of incidence
- Lift coefficient for aerodynamic models in supersonic flow
- Shock waves and expansion patterns around a two-dimensional model in supersonic flow conditions (when used with the optional Schlieren Apparatus).

Recommended Ancillaries

- Schlieren Apparatus AF300a
- Air Compressor Receiver and Dryer AF300b

Essential Ancillaries

- Versatile Data Acquisition System - Frame-mounted version (VDAS-F)

Essential Services

Electrical supply:

50 Hz to 60 Hz 90 VAC to 250 VAC for the AF300 Instrument Frame. 'See specifications' for AF300b.

Air supply:

Greater than 7 bar and $0.5 \text{ kg}\cdot\text{s}^{-1}$ of clean, dry air, or use the optional AF300b.

The air supply must not contain more than 0.2 kg of water for 450 kg of air.

Floor space needed:

4 m x 2 m for the wind tunnel and optional Schlieren Apparatus (AF300a)

5 m x 3 m for the wind tunnel with the optional Compressed Air Supply (AF300b)

Operating Conditions

Operating environment:

Laboratory environment that allows for high sound levels

Storage temperature range:

-25°C to +55°C (when packed for transport)

Operating temperature range:

+5°C to +40°C

Operating relative humidity range:

80% at temperatures < 31°C decreasing linearly to 50% at 40°C

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecequipment.com • **W** www.tecequipment.com
- An ISO 9001 certified company
- VDAS is a registered trademark of TecEquipment Ltd



TecEquipment

AF300

Intermittent Supersonic Wind Tunnel

Sound Levels

This equipment emits sound levels greater than 100 dB(A). You must wear ear defenders when you use it or work near to it.

Specifications – AF300

Dimensions and weight (wind tunnel):

Nett: 2000 mm high x 3031 mm long x 805 mm wide and 236 kg

Packed: Approximately 5.9 m³ and 300 kg

Dimensions and weight (instrument frame without instruments):

Nett: 1260 mm long x 840 mm high x 510 mm wide and 22 kg

Packed: Approximately 0.78 m³ and 30 kg

Nominal nett dimensions (AF300b):

2120 mm high x 4500 mm x 1000 mm

Working section:

Nominal 100 mm x 25 mm

Liners:

Subsonic, Mach 1.4 and Mach 1.8

Models:

- 5-degree single wedge
- 7-degree double wedge
- 10-degree double wedge
- 10-degree double wedge with two pressure tappings.

Instruments:

- Angle encoder input board for VDAS-F
- Angle encoder
- 32-way pressure display
- Pressure mimic module
- Delivery pressure – mechanical gauge and electronic transducer
- Supply pressure – mechanical gauge

Specifications – AF300b

- Three air receivers/reservoirs – one receiver has its own inlet and outlet isolation valves to give a choice of operating time and pressure recharge rate
- Main outlet isolation valve
- Maximum pressure approximately 14 bar
- Compressor – rated at a nominal 15 kW total electrical power, needing 32 A at 400-440 VAC three phase with earth. Voltage to be specified on order.
- Air dryer and filter – rated at less than 25 W and normally single phase. Voltage to be specified on order.

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecequipment.com • **W** www.tecequipment.com
- An ISO 9001 certified company
- VDAS is a registered trademark of TecEquipment Ltd

tradition.

innovation.

integration.

infoWERK is a leading expert in the development of eLearning courseware, learning system solutions, teaching and AV equipment.

Furthermore infoWERK is the representative and system integrator of "TecQuipment".

TecQuipment is one of the global leaders in technical teaching equipment for engineering. If you are interested in one of TecQuipment's products feel free to contact us at:



infoWERK Medien & Technik GmbH

Martinsbühel 6 / A-6170 Zirl / Austria

Phone: +43 (0) 5238 52099-0 / Fax: +43 (0) 5238 52099-40

E-Mail: info@infowerk.at / Website: infowerk.at

Otto-Dürr-Straße 25

D-70435 Stuttgart, Zuffenhausen/ Germany

Phone: +49 (0) 711 342471-0 / Fax: +49 (0) 711 342471-11

E-Mail: info@de.infowerk.at / Website: infowerk.at