




Aerodynamics

AF13

Round Turbulent Jet

Allows students to investigate a jet of air as it emerges from the end of a tube

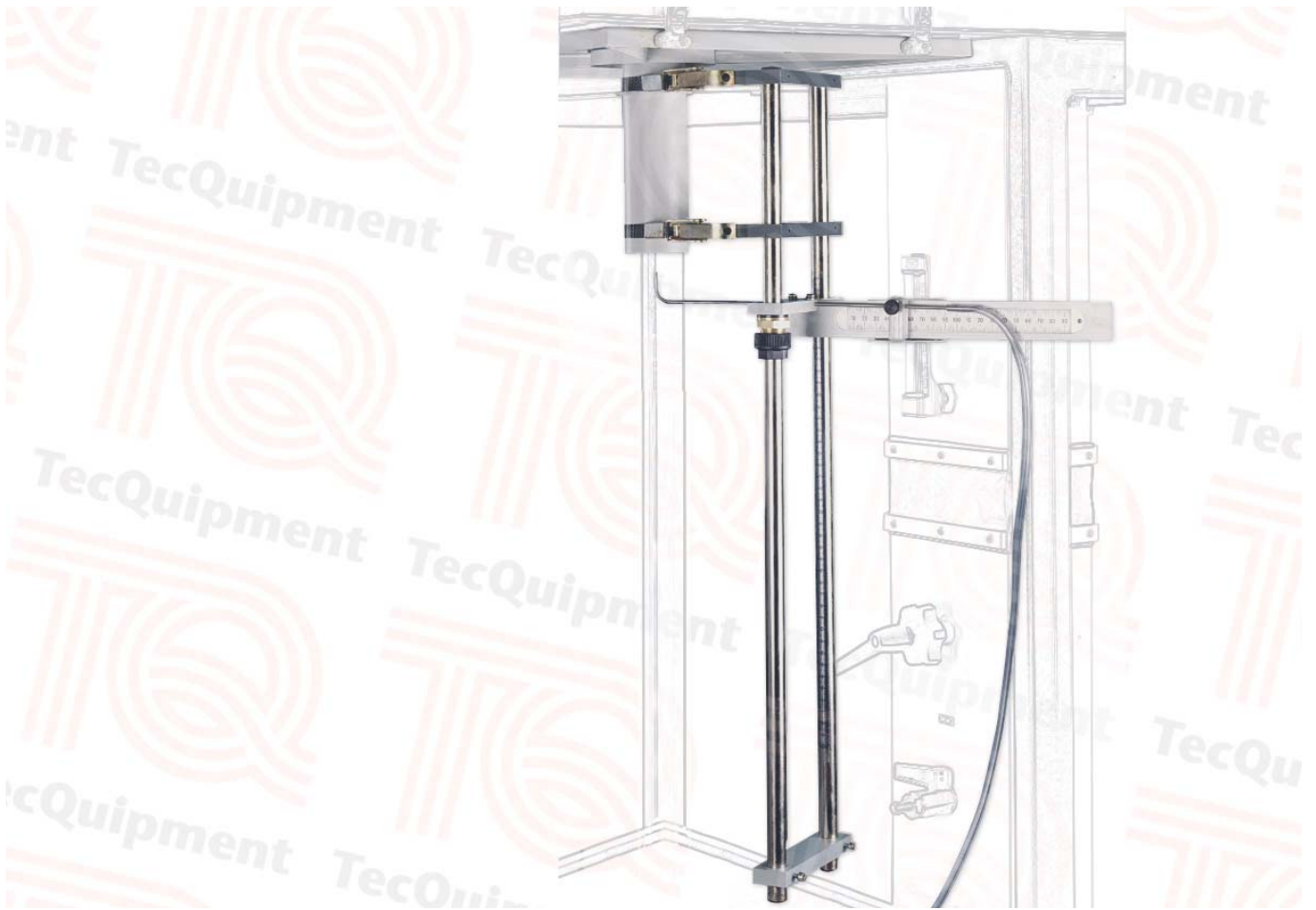

- One of a series of eight experiment modules that fits to the Modular Air Flow Bench (AF10)
- Allows a number of tests on the velocity of a submerged jet emerging from the end of a tube
- The tube has a carefully designed inlet for best results
- Toggle clamp connections to the Modular Air Flow Bench plenum chamber for quick and easy fitment
- Quick-release coupling fitted to the Pitot tube to allow rapid and reliable connection to the AF10a Manometer

AF13

Round Turbulent Jet

Description

This module consists of a tube with a specially designed rounded entry. The tube is mounted in a stiff plate with the rounded entry on one side and the exit on the other.

To set the experiment up the contraction is unclipped from the Air Flow Bench (AF10) and set aside. The whole plate is then mounted onto the plenum chamber directly so that air enters the rounded entry of the tube and leaves the end in a jet. The total pressure in the jet is measured by a Pitot tube held by a traversing mechanism. The mechanism allows the Pitot tube to move radially across the jet and axially along it. The Pitot tube is connected to the AF10a manometer (ancillary) via a flexible tube fitted with a quick-release coupling.

Standard Features

- Supplied with a comprehensive User Guide
- Five-year warranty
- Manufactured in accordance with the latest European Union directives

Essential Base Unit

- Modular Air Flow Bench (AF10)

Essential Ancillaries

- Multitube Manometer (AF10a)

Experiments

- Decay of the centre-line velocity.
- Velocity profile at various distances along the jet and the development of the spread of the jet.
- Analysis of the velocity profiles to show how the mass flux in the jet increases, the kinetic energy flux decreases and the momentum flux remains constant along the jet length.

Specifications

Packed dimensions and weight:
0.08 m³; 15 kg

Operating Conditions

Operating environment:
Laboratory

Storage temperature range:
–25°C to +55°C (when packed for transport)

Operating temperature range:
+5°C to +40°C

Operating relative humidity range:
80% at temperatures < 31°C decreasing linearly to 50% at 40°C

- TecEquipment Ltd, Bonsall Street, Long Eaton, Nottingham NG10 2AN, UK
- **T** +44 115 972 2611 • **F** +44 115 973 1520 • **E** info@tecequipment.com • **W** www.tecequipment.com
- An ISO 9001 certified company

tradition.

innovation.

integration.

infoWERK is a leading expert in the development of eLearning courseware, learning system solutions, teaching and AV equipment.

Furthermore infoWERK is the representative and system integrator of "TecQuipment".

TecQuipment is one of the global leaders in technical teaching equipment for engineering. If you are interested in one of TecQuipment's products feel free to contact us at:



infoWERK Medien & Technik GmbH

Martinsbühel 6 / A-6170 Zirl / Austria

Phone: +43 (0) 5238 52099-0 / Fax: +43 (0) 5238 52099-40

E-Mail: info@infowerk.at / Website: infowerk.at

Otto-Dürr-Straße 25

D-70435 Stuttgart, Zuffenhausen/ Germany

Phone: +49 (0) 711 342471-0 / Fax: +49 (0) 711 342471-11

E-Mail: info@de.infowerk.at / Website: infowerk.at